Prism Spectrometer

Rahmanyaz Annyyev, Hikmat Gulaliyev 02 March 2024

Abstract

In this experiment, we utilize a prism spectrometer—an optical device that separates light into its constituent frequencies—to measure the refractive index of a prism. A collimator with a slit is used to produce a parallel beam of light emitted by a mercury lamp, which is then passed through the prism. The light is then refracted and dispersed into its constituent frequencies. The angle of deviation of the light is measured and used to calculate the refractive index of the prism using a special relation. A graph of the index of refraction versus wavelength and the inverse square of the wavelength is plotted, and the constants of the Cauchy equation are obtained. The results are in line with the theoretical values, and the experiment is successful.

1 Introduction

The visible spectrum is the portion of electromagnetic radiation that is visible to the human eye. It is composed of light with wavelengths between 380 and 780 nm. Each wavelength, or frequency of light is perceived as a different color by the human eye [1]. Spectrum analysis is the study of measure and interpretation of electromagnetic spectra. It is used in various fields such as physics, chemistry, and astronomy. One of the most important tools in spectrum analysis is the prism spectrometer.

A prism spectrometer is an optical device, and it is composed of a collimator, a prism, and a telescope. The collimator is a device that produces a parallel beam of light. It is composed of an adjustable slit on one end and a converging lens on the other. The slit is used to control the width of the entering light, and the lens is used to focus the light into a parallel beam. As the beam exits the collimator, it impinges on the prism. The prism rests on a rotatable table, and it is used to refract and disperse the light into its constituent colors. The telescope is used to observe the dispersed light, and it is composed of a lens and an eyepiece. The set of parallel beams of different wavelengths form images of the slit. The telescope can be rotated with the aid of the reticle to position the images for each spectral line. The instrument that is used throughout the experiment to measure the angles is called a vernier scale. Note that a mercury lamp was used as a source of light.

The experiment is divided into two parts. Before proceeding, we follow two preliminary steps. First, we adjust the telescope on a distant object to ensure that it is properly focused. Next, we calibrate the collimator to ensure that the light is parallel.

In the first part, we measure the apex angle of the prism—the angle between the two faces of the prism. The prism is placed on the table such that the angle is towards the collimator. Next, we rotate the telescope and find two positions where the light is *reflected* from the prism. The image of the slit is observed through the telescope, and the telescope is further adjusted until the image falls on the reticle. The positions are then read from the vernier scale. The apex angle is given by the relation

$$\alpha = \frac{\text{Position1} + \text{Position2}}{2}, \tag{1}$$

where $Position_1$ and $Position_2$ are the two angles. Several measurements are taken, and the average is calculated using the formula

$$\bar{\alpha} = \frac{1}{N} \sum_{i=1}^{N} \alpha_i, \tag{2}$$

where N is the number of measurements.

In the second part, we use the spectrometer to measure the refractive index of the prism. It is a measure of how much the light is bent as it passes through the prism. Examine Figure 1. As you can see, δ is the angle of deviation of the light—the angle between the incident and emergent beams of light, and α is the angle of the prism. To find the refractive index of the prism, we measure the minimum deviation angle of the light.

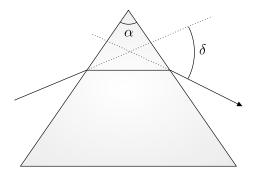


Figure 1: Refraction of light through a prism.

The minimum angle of deviation occurs when the angle of incidence is equal to the angle of emergence. Hence, in this case, the refractive index of the prism is given by the relation

$$n = \frac{\sin((\delta_{min} + \alpha)/2)}{\sin(\alpha/2)},$$
(3)

which is derived from Snell's law [3],

$$n_1 \sin(\theta_1) = n_2 \sin(\theta_2),\tag{4}$$

where n_1 and n_2 are the refractive indices of the two media, and θ_1 and θ_2 are the angles of incidence and refraction, respectively. The methodology for finding the minimum deviation angle is as follows. We need to determine two positions. First, we remove the prism from the table and record the position of the telescope such that it is parallel to the collimator and the image of the slit is adjusted to fall on the reticle. This is the first position. Next, we place the prism back on the table such that one of the faces makes an angle of 45° with the beam coming from the collimator. Next, we rotate the telescope and find the position where the light is refracted from the prism. The spectral lines in several colors are observed through the telescope, and the latter is further adjusted until each line falls on the reticle. Then, the telescope is rotated further until a point where the direction of movement of the spectral lines reverses in direction. This is the second position. The deviation angle is given by the relation

$$\delta_{min} = 360^{\circ} - \text{Position1} + \text{Position2}.$$
 (5)

A dispersion curve is a graph of the refractive index versus wavelength. It can be represented using the Cauchy equation [2],

$$n(\lambda) = A + \frac{B}{\lambda^2} + \frac{C}{\lambda^4} + \cdots, \tag{6}$$

where A, B, and C are constant pertaining to a specific material. The dispersion curve is used to determine the refractive index of a material at a specific wavelength.

2 Data & Results

The results of the first part of the experiment are shown in Table 2. The average angle of the prism is given by $\bar{\alpha} = 59.87^{\circ}$.

Position #1	Position #2	α_i	
52°46′	293°01′	59.86°	
52°52′	293°05′	59.88°	
52°48′	293°03′	59.86°	
		$\bar{\alpha} = 59.87^{\circ}$	

Table 1: Data for the prism angle, α .

The results of the second part of the experiment are shown in Table 2. The wavelengths of the spectral lines were taken from the table given in the laboratory manual.

λ (nm), Color	Position #1	Position #2	δ_{min}	n
404.6, Violet	351°48′	27°40′	40.19	1.536
435.8, Blue	351°48′	28°41′	39.18	1.524
546.1, Green	351°48′	29°28′	38.31	1.514
579.1, Orange	351°48′	29°49′	38.02	1.511

Table 2: Data for the refractive index versus wavelength.

After evaluating the refractive index for each wavelength, two graphs were plotted. The results are shown in Figures 2 and 3. The first graph shows the relationship between the refractive index of the prism and the wavelength of the light. The second graph shows the relationship between the refractive index and the inverse square of the wavelength, and the line is the least squares fit of the data. The values of the constants A and B are obtained: A = 1.58 and B = 0.00734.

$$n(\lambda) = 1.58 + \frac{0.00734}{\lambda^2}. (7)$$

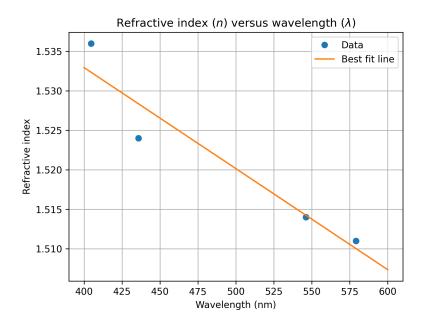


Figure 2: Refraction index versus wavelength.

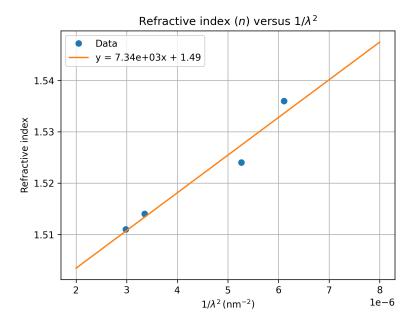


Figure 3: Refraction index versus inverse square of wavelength.

3 Discussion & Conclusion

Due to the analog nature of this experiment, human error is inevitable. In particular, manually reading the angles using a vernier scale introduces some errors. Moreover, manually adjusting the telescope to find the minimum deviation angle is another source of error.

Before and during the experiment, certain assumptions were made, including the assumption that the light source was monochromatic, the prism was perfectly shaped and uniform, the refractive index of air was 1, and the telescope was perfectly calibrated. However, these assumptions may not hold in real-life scenarios and could introduce errors in the results.

Despite these limitations, the experiment aligns with the theoretical values in every aspect except for the number of spectral lines. This discrepancy might be explained by the fact that some of these lines are too close together for the human eye to discern.

In conclusion, the experiment was successful, with results largely aligning with the theoretical values. Introducing optoelectronics in further studies could minimize human error and enhance the accuracy of the results.

References

- [1] Robert T. Marcus. *The Measurement of Color*. Vol. 1. AZimuth. North-Holland, 1998, pp. 31–96. ISBN: 978-0-444-82825-6. DOI: https://doi.org/10.1016/S1387-6783(98)80005-6.
- [2] Frank L. Pedrotti and Leno S. Pedrotti. *Introduction to Optics*. Pearson Education, 2013, pp. 63–64. ISBN: 9780131499331.
- [3] Eugene Hecht. *Optics*. 5th ed. Pearson Education, 2017, pp. 108–110. ISBN: 1-292-09693-4.